

AURA REVIVE

Deep Muscle Pain Relief Device



 powered by
auramedical

USER MANUAL

Model No. AR40

CAUTION

Do not immerse in water. Keep away from liquids or heat sources.
Keep ventilation ports free from dust and debris.
Do not remove screws or attempt to disassemble.
Do not operate continuously for more than one hour.
Allow the device to rest for 30 minutes before reusing.
Unplug the charger from the device after charging or prior to use.
For adult use only. Do not use if injured.
Consult your doctor before using this product.

SAFETY INSTRUCTIONS

WARNING: TO REDUCE RISKS OF ELECTRIC SHOCK, FIRE AND PERSONAL INJURY, THIS PRODUCT MUST BE USED IN ACCORDANCE WITH THE FOLLOWING INSTRUCTIONS!

Use only on a dry clean surface of the body, over clothing, by lightly pressuring and moving across for approximately 60 seconds per region.

Use the Aura Revive only on the body's soft tissue, as desired, without causing pain or discomfort.

Do not use on your head or any hard or bony portion of your body.

Use only the accessory heads that provide the best desired effect.

Bruising can occur regardless of the control setting or pressure applied.

Check treated areas frequently and immediately stop at the first sign of pain or discomfort.

Keep fingers, hair or any other body part away from the shaft and backside of the device as pinching may occur.

Do not place any objects into the ventilation holes of the device. Do not immerse in water or allow water to enter the ventilation holes.

Do not drop or misuse the device. Only charge with the supplied charger.

Carefully examine the device and battery before each use. Do not tamper with or alter the device in any way. Never leave the device unattended while charging or operating.

Use a slightly damp towel to clean the device and dry with a soft cloth.

WARNING

PLEASE DO NOT USE THE DEVICE OR ANY VIBRATION TOOL WITHOUT OBTAINING MEDICAL APPROVAL!

Use cautiously and consult a medical professional before using with pregnancy, diabetes with complications such as neuropathy or retinal damage, pace-makers, recent surgery, epilepsy or migraines.

Use caution with herniated disks, spondylosis or similar conditions, recent joint replacements or LUD's, metal pins, plates or any concerns about your health. Frail individuals and children should be accompanied by an adult when using any vibration device.

Ongoing research is being done regarding the effects of vibration exercises for specific medical disorders. This will likely result in the shortening of the list of contraindications mentioned above.

Practical experience has shown that integration of vibration exercises into a treatment plan is advisable in a number of cases. This must be done on the advice of and in the presence of a doctor, specialist or physiotherapist.

WARRANTY

IF WITHIN ONE YEAR FROM THE DATE OF PURCHASE THIS PRODUCT FAILS DUE TO A DEFECT IN MATERIAL OR WORKMANSHIP, WE WILL REPAIR OR REPLACE THE PRODUCT OR NECESSARY COMPONENTS, FREE OF CHARGE.

THE WARRANTY EXCLUDES:

Damage caused by accident, abuse, mishandling or transport.

Units subjected to unauthorised repair.

Units not used in accordance with the instructions.

Damage exceeding the cost of the product.

Deterioration of the delivered product resulting from abnormal storage and/or safeguarding conditions by the client.

Failure to provide a dated proof of purchase.

**Manufactured in China for Aura Medical LLC
888.433.6916 | support@aura-medical.com**

THANK YOU FOR PURCHASING THE AURA REVIVE!

The Aura Medical Team develops and manufactures high end medical and health focused devices designed to effect positive change in people's lives. We always strive for innovation. We work hard to be a step ahead of the current market leaders in terms of design, technology, consumer health and satisfaction.

SPECIFICATIONS



Vibration Range: 12mm
Speed Range: 1100 - 3500rpm
Decibel Range: 40 - 60db
Rated Voltage: 24VDC
Rated Power: 18W
Speed: 5 speeds adjustable
Charging Voltage: 24VDC
Input Voltage Frequency: 50 - 60Hz
Rated Output Current: 0.5A
Plug: US/UK/AU etc.
Warranty: 12 months
Certification: CE/FCC/RoHS
Dimensions: 160 x 213 x 60mm
Weight: 1.1kg
Accessories: 8pcs.
Chargin Time: 2 - 3hrs.
Operation: Approx. 3hrs.
Battery Capacity: 2600mAh
Model no: AR40

USING THE DEVICE

Before initial use we recommend charging the device for about 3 hours. To recharge, please connect the unit with the supplied AC power cable to the cable socket on the unit.

We recommend recharging the device when the battery power level display is showing only 1 bar lighting up or when the LED indicator lights up red.



RED:
Power Low

Please recharge the device



GREEN:
Power Ready

FLASHING GREEN:
Battery is recharging



- VIBRATION SPEED DISPLAY**
- HEATING FUNCTION DISPLAY**
- STANDARD PRESSURE MODE**
- ASCENDING PRESSURE MODE**
- FLUCTUATING PRESSURE MODE**
- POWER BUTTON**

Align and insert the battery into the base of the unit. Turn on the device by moving the switch at the bottom of the battery. When the switch is in the ON position the device can be switched to the first intensity level by pressing the POWER button on the back of the device. To increase the intensity continue pressing the POWER button as needed.

The device has a unique heat feature for deep muscle relief. Please note that the heating feature is only available with the HEAT HEAD accessory. Insert the head and ensure a tight fit with the unit. Power on the device and press the STANDARD PRESSURE mode TWICE or until the HEATING FUNCTION light turns on.

MAIN SWITCH AND CHARGING SWITCH

ACCESSORIES



Heating Head
can promote
blood circulation
more effectively



Fork Head
suitable for
both sides of
the spine




Lotus Head
suitable for
the soles of
your feet



Round Head
suitable for
the shoulders
and neck



Shovel Head
suitable for
the spine

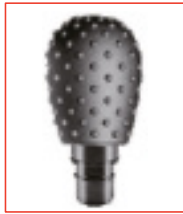

Mace head
suitable for large
muscles such as
hips and legs



Bullet Head
suitable for
local pain points
and small muscle groups



Five-Claw Head
for muscles
bulging such
as calf muscles

AURA REVIVE